

hipocross 2.0

Brace for hip and coxalgia



ORTHOSERVICE
RO+TEN

Take care feel better

Several applications, several benefits

APPLICATIONS

BENEFITS

post-arthroscopy care	▶ oedema prevention, cryotherapy and movement control
conservative treatment of hip dislocations	▶ movement control (flexion and extension and especially postoperative abduction/ adduction)
post-surgical care (also for replacements or revision surgery)	
arthrosis of the hip joint	▶ pain relief
coxalgia or inoperable joint misalignment	

INNOVATION 2

Only 3 fits for any request

The versatility of the brace reduces the need to have to keep stocks in stock, knowing nonetheless that every client's request will be satisfied.

Size	S	M	L
Pelvis circ. cm	72/104	104/137	137/170
Proximal thigh circ. cm	40/58	54/74	68/90

Already a winner, but even better from today: 3 absolute innovations

TELESCOPIC AND ARTICULATED ROD

INNOVATION 1

Allows numerous regulations:

- **Height**, by sliding the two segments to the desired measurement.
- **Adduction**, thanks to the joint present at trochanteric height
- **Flexion-extension**, through a special graduated joint (from 0° to 120°).



INNOVATION 3



VELCRO® FASTENING SYSTEM FOR THE STRAPS ON BOTH THE THIGH AND THE HIP BONE

The Velcro® allows quick and millimeter precision fastening.



ANTI-EDEMA CONTAINMENT FABRIC AND CRYOTHERAPY GEL


The compression exerted by the combination of the special containment fabric and the gel (that has to be inserted in the pocket in the area of the hip), help the reabsorption of the edema. The gel is cooled in the freezer.

Attention to detail
makes the difference



RIGID BAR

A rigid bar that locks the joint in flexion and neutral abduction during the initial days following a dislocation.

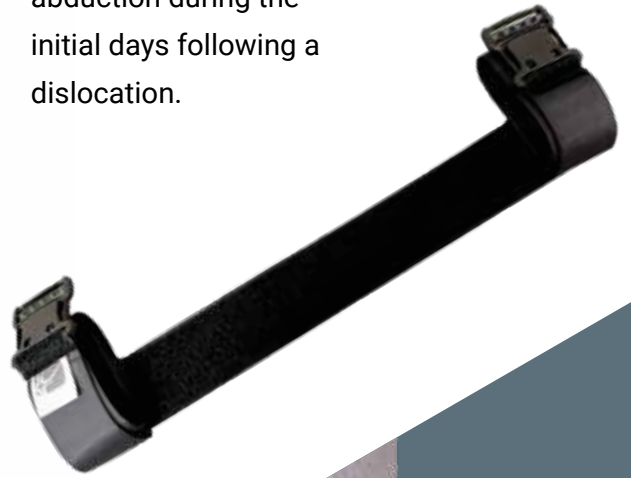

X-SHAPED TIE RODS ON THE SIDE FOR EASY ADDUCTION

The 2 robust rods that cross on the side allow the dynamic adjustment of the adduction.



hipocross 2.0

Brace for hip and coxalgia

REF.91006N - REF.91007N

Indications

REF.91006N:

- Post-arthroscopy care.
- Post-surgical care (compatible with the degree of stabilization required).
- Conservative treatment of hip dislocations.
- Arthrosis of the hip joint.
- Coxalgia or inoperable joint misalignment.



REF.91007N

Adjustable (telescopic) rod for Hipocross



SINGLE PATIENT
MULTIPLE USE



PATENT
EP3202376



LOW
PROFILE



BREATHABLE



3D FABRIC



ANATOMICAL
FIT



TELESCOPIC

Item	REF.91006N (Hipocross 2.0)		
Size	S	M	L
Pelvis circ. cm	72/104	104/137	137/170
Proximal thigh circ. cm	40/58	54/74	68/90
Colour	black		

specify right or left

Item	REF.91007N (Adjustable telescopic rod for Hipocross 2.0)
Size	universal
Colour	black

specify right or left



orthoservice.com



Headquarter: ORTHOSERVICE AG
Via Milano 7 · CH-6830 Chiasso (TI) · Switzerland
Tel. 0041 (0) 91 822 00 88 · Fax 0041 (0) 91 822 00 89
info@orthoservice.com · www.orthoservice.com

Niederlassung Deutschland: Orthoservice Deutschland GmbH
Flugstraße 8 · D-76532 Baden-Baden · Deutschland
Tel. 0049 (0) 7221 991 39 11 · Fax 0049 (0) 7221 991 39 13
info@orthoservice.de · www.orthoservice.de

Siedziba w Polsce: Orthoservice Polska Sp. z o. o.
ul. Warszawska 416a · 42-209 Częstochowa · Polska
Tel. +48 (0) 34 340 13 10
info@orthoservice.pl · www.orthoservice.pl

Sede italiana: RO+TEN S.r.l.

Sede legale: Via Marco De Marchi, 7 · I-20121 Milano (MI) · Italia

Sede operativa e amministrativa:
Via Comasina, 111 · I-20843 Verano Brianza (MB) · Italia
Tel. 0039 039 601 40 94 · Fax 0039 039 601 42 34
info@roplusten.com · www.roplusten.com

Società soggetta a Direzione e Coordinamento (art. 2497bis CC):
Orthoservice AG (CH) · 6830 Chiasso (TI) · Switzerland



**ORTHOSERVICE
RO+TEN**

Take care feel better