

lumbopull 70/72

**Intuition and engineering
against backache**

NEW



**ORTHOSERVICE
RO+TEN**

Take care feel better



With Lumbopull you can...

Lumbopull achieves
all these goals because:


- 1 It **increases** intra-abdominal pressure
- 2 It **promotes** spinal decompression
- 3 It **creates** a natural barrier to flexion,
extension, and rotation
- 4 It **preserves** the movements necessary
for performing everyday activities

To contain is to support, that's the secret...

The **lumbar section of the spine** is a critical point, often the site of pain of varying intensity.

Whether it's a **transient low back pain** or a **post-operative recovery**, the elements that best promote healing are:

- 1 Reduction of the load on the vertebrae
- 2 Correct posture
- 3 Limitation of movement

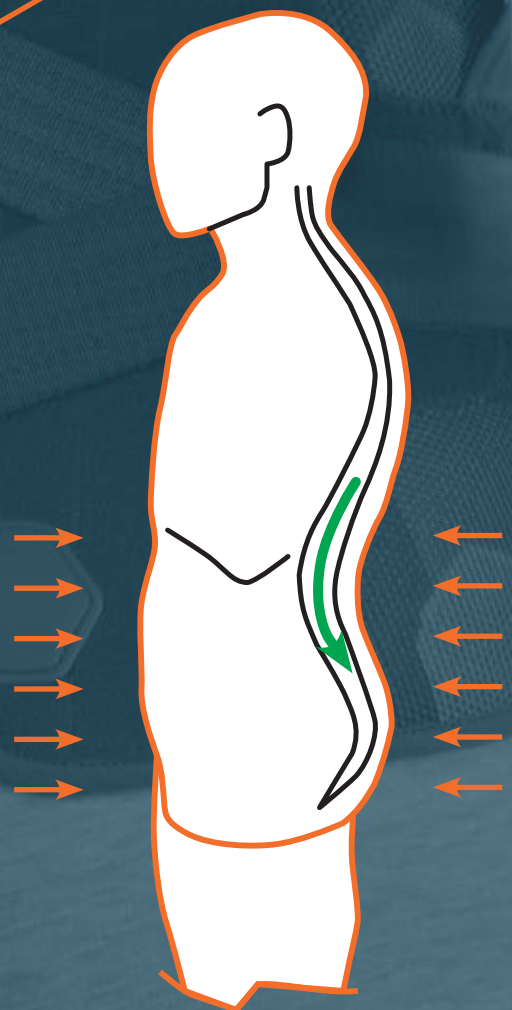

WHIT LUMBOPULL

lumbopull 70/72

**Quadruple closure strength
for maximum support**



Exclusive strap a

Maximum traction with
to make fastening easier even f



EXCELLENT FIT AND ADAPTABILITY

with attention
to style

ANATOMICAL DESIGN

perfectly follows
the body's contours
and makes the brace
easy to wear



ONE FABRIC, TWO FUNCTIONS:

- **containment** through the external technical "blocked" fabric
- **breathability** thanks to the internal three-dimensional mesh

and pulley system

minimal effort required
for weakened or elderly patients



Rear strap pass-through to ensure adaptation to waist-pelvis proportions



SMALL LUMBAR PAD

made of polymer material to ensure **stability** and **rigidity**, thermoformable



REAR ELASTIC INSERT

made of soft fabric, allows Velcro® application of various options

Effective with an exceptionally low profile

Dorsopull 74: the brace becomes a shoulder harness in just a few steps

In pathological conditions involving the **thoracic spine**, the brace can become a shoulder harness simply by adding the optional **Dorsopull** back panel



- **ADJUSTABLE SHOULDER STRAPS**

with elastic return at the back to maximize the applied force

- **MADE OF PADDED FABRIC**

- **SELF-MOLDING STAYS**

for increased internal strength



- **RIGID FRAME**

made of aluminum alloy, **modal**e and **radiolucent**. Can be removed when necessary



lumbopull 70 | REF.90070

Low orthosis with easy-pull tensioning straps for stabilising and supporting the lumbar spine, with rigid thermoformable pad

Code	REF.90070				
Size	S	M	L	XL	XXL
Pelvic Circumference (cm)	70/85	85/100	100/115	115/130	130/145
Front Height (cm)	16	16	16	16	16
Back Height (cm)	26	26	26	26	26

Color: grey

Option: REF. 70020 Large lumbar pad 24x31



REF. 90070



REF. 90072

lumbopull 72 | REF.90072

High orthosis with easy-pull tensioning straps for stabilising and supporting the lumbar spine, with rigid thermoformable pad

Code	REF.90072				
Size	S	M	L	XL	XXL
Pelvic Circumference (cm)	70/85	85/100	100/115	115/130	130/145
Front Height (cm)	20	20	20	20	20
Back Height (cm)	31	31	31	31	31

Color: grey

Option: REF. 70020 Large lumbar pad 24x31



REF. 90074

Option: REF.70020
Large lumbar pad**dorsopull 74 | REF. 90074**

Padded shoulder back splint for the orthopedic Lumbopull orthosis

Code	REF.90074			
Size	XS	S	M	L
Distance T1-L5 (cm)	35/41	41/47	47/54	54/61

Color: grey

Indications**LUMBOPULL 70/72**

- Lumbago
- Sciatic lumbalgia and lumbar cruralgia
- Spondylarthrosis and lumbar discopathies
- Slight traumas to the lumbosacral column and paravertebral lumbar contractures
- Post-operative recovery
- After-effects of stable compression fractures of lumbar vertebrae
- Treatment of fractures to transverse processes
- Only Lumbopull 72: Degenerative lumbar scoliosis
- Only Lumbopull 72: Spondylolysis with or without vertebral listhesis

DORSOPULL 74

- Back pain due to arthrosis or kyphosis
- After-effects of fractures of dorsal vertebrae and/or upper lumbar vertebrae
- After-Effects of stable compression fractures of the dorsal spine due to osteoporosis or metastasis

Contraindications

- Currently none known



www.orthoservice.com



Headquarter: ORTHOSERVICE AG
Via Milano 7 · CH-6830 Chiasso (TI) · Switzerland
Tel. 0041 (0) 91 822 00 88 · Fax 0041 (0) 91 822 00 89
info@orthoservice.com · www.orthoservice.com



ORTHOSERVICE DEUTSCHLAND GmbH
Flugstraße 8 · D-76532 Baden-Baden · Deutschland
Tel. +49 (0) 7221 40 533 0 · Fax 0049 (0) 7221 40 533 20
info@orthoservice.de · www.orthoservice.de

ORTHOSERVICE POLSKA Sp. z o.o.
ul. Warszawska 416a · 42-209 Częstochowa · Polska
Tel. +48 (0) 34 340 13 10 · (NIP): 9492246785
info@orthoservice.pl · www.orthoservice.pl

OSUK Ltd
Units 47&50 North East BIC
Wearfield
Sunderland Enterprise Park East
Sunderland SR5 2TA
Tel +44 (0) 191 51 66 220
info@orthoservice.uk · www.orthoservice.uk

RO+TEN S.r.l.
Sede legale: Via Marco De Marchi 7 · I-20121 Milano (MI) · Italia
Sede operativa e amministrativa:
Via Comasina 111 · I-20843 Verano Brianza (MB) · Italia
Tel. 0039 0362 18 51 200
info@roplusten.com · www.roplusten.com



PATENTED



EASY PULLING



BREATHABLE



ANATOMICAL FIT



3D FABRIC



SINGLE PATIENT MULTIPLE USE



**ORTHOSERVICE
RO+TEN**

Take care feel better